


FileCare® EDTA **chelating agent**

Better results for root canal treatments

→ Why FileCare® EDTA ?

- **removes smear layer**
- **improved tubular penetration**
of endodontic medication and filling material,
as dentinal tubules are opened
- **helps prevent files from blocking,**
- **reduces instrument fracture**
reduces friction on the instrument, because
FileCare® EDTA works as a lubricant and
pre-softens dentin

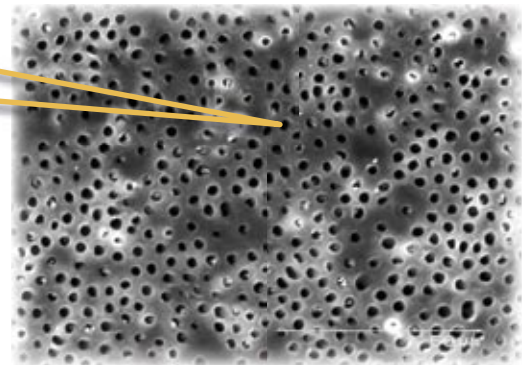

coronal portion of canal wall after preparation
original magnification 1000x



- **easier penetration of the instrument,**
as FileCare® EDTA decalcifies dentin by binding the
calcium in an electro-negative complex, thus
overcoming blockages within the root canal.
- **pulp tissue and dentinal debris are
actively flushed out of the canal,**
due to the foaming effect of FileCare® EDTA
in contact with NaOCl
- **anti-bacterial effect**
from EDTA and urea peroxide

More than just a lubricant



FileCare® EDTA

FileCare® EDTA is a chelating agent which reduces stress on the file and allows efficient cleaning and easier preparation

simple to use



Simply put some FileCare® EDTA from a dappen dish or similar...



...onto an endodontic instrument (e.g. FlexMaster® or hand-file)...



...and introduce it into the canal.

FileCare® EDTA contains:

15% EDTA

10% urea peroxide gel



Ref. 1010 000 002 pack with **two** 3ml syringes

Ref. 1010 000 005 pack with **five** 3ml syringes



VW000130 Rev. 0/3.12.2002

VDW
ENDODONTIC SYNERGY

VDW GmbH
P.O. Box 830954 • 81709 Munich, Germany
Tel. +49 89 62734-0 • Fax +49 89 62734-195
info@vdw-dental.com • www.vdw-dental.com

Endo Easy Efficient®